

Brannan instruction sheet

BODY TEMPERATURE INDICATOR


Product ref: 11/466/9

REF F3



MD

CE 0051

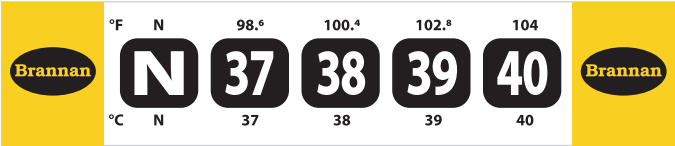

Instructions for use ():

Body temperature indicator with encapsulated liquid crystals (). Place the thermometer for 15 seconds on a dry forehead, holding it by both ends. If the letter N is visible there is no fever. If a number appears, the temperature corresponds to the visible value. If two numbers are visible, the temperature is a value in between. If temperature rises, consult a doctor. Use the indicator at room temperature, with adequate brightness, under adult supervision for children under 13. Clean the indicator by rubbing it gently with soapy water for at least 30 seconds or with alcohol free skin disinfectants (do not soak). External use on healthy skin. Do not expose to sunbeams (), keep away from heat sources, alcohol and solvents. (Issue of June 17, 2020).

 **Made for Brannan by:**
I.P.S. S.r.l. - Via Civesio, 6 - 20097 San Donato Mil.se (Mi) - Italy +39 025279641
www.ips-srl.it

Explanation of symbols:

CE 0051 CE mark and identification of notified body

 Manufacturer's name and date of production of the device

 Product lot number

 Device model

 Consult instructions for use

 Do not expose to sunlight or UV sources

 Do not dispose of in the environment. Dispose of consciously.

 Polycoupled material. Mixed waste collection

°C Degrees Centigrade

°F Degrees Fahrenheit

 Medical Device

