

Product instruction sheet

Multi-Function Timer (28/219/0)

Thank you for choosing this high quality, multi-function digital timer.

- Please read this instruction manual carefully to familiarise yourself with your new timer.
- Carefully following the instructions in this manual will prevent damage to your instrument and loss of your statutory rights arising from defects due to incorrect use.
- Please keep this instruction manual for future reference.

Functions:


- Clock • Stopwatch • Timer • Vibration, audible and visual alarms • Clip, stand and magnet fixtures

Timer Description:

- (1) Light signal
- (2) LCD display
- (3) MIN and SEC keys
- (4) START/STOP key



- (5) Light signal ON/OFF switch
- (6) Audible signal ON/OFF switch
- (7) Power ON/OFF switch
- (8) Vibration signal ON/OFF switch
- (9) Magnet
- (10) Fastening clip

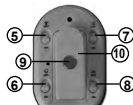

Preparing your Timer

1. Remove all packaging and dispose of appropriately.
 2. Open the battery compartment on the reverse of the timer by pushing the cover with the fastening clip downwards whilst applying small pressure.
 3. Remove the security tab to activate battery.
 4. Close the battery compartment.
- The timer is now ready for use.

Setting the Clock

1. Switch on the timer using the power on/off switch. The display will show "12:00".
2. Hold [MIN] for 2 seconds, the hour time will flash. Set the desired hour by pushing [SEC].
3. Press [MIN] again (do not hold down) then [SEC] to set the minutes.
4. Press [MIN] again to confirm time.

Setting the Stop Watch

1. Hold [MIN] for 4 seconds to switch from clock to timer display. Initially minutes/seconds (M/S) are displayed with a start time of "00:00".
 2. Press [START/STOP] to activate the stop watch function. This can be interrupted and continued by pressing the [START/STOP] again. The maximum measurable time is 99 minutes and 59 seconds.
 3. To set the stop watch for hours/minutes (H/M), push [MIN]. The time is measured as described above. The maximum measurable time is 99 hours and 59 minutes.
- Press [MIN] again to switch from timer to clock.

Setting the Timer

Note: The time lapsed can be indicated with a combination of up to three alerts: an audible signal, a vibration alert and/or a light signal by setting the corresponding switch on the reverse side of the timer to position "ON".

1. Hold [MIN] for 4 seconds to switch from clock to timer display. Initially minutes/seconds (M/S) are displayed with a start time of "00:00". Set the desired minutes using the [MIN] button and the desired seconds using [SEC] button.
2. Press [START/STOP] to activate the timer. This can be interrupted and continued by pressing the [START/STOP] again. The maximum measurable time is 99 minutes and 59 seconds.
3. When the required time has elapsed the timer blinks, beeps or vibrates, according to the chosen signal.
4. Stop the alarm by pushing any button.
5. To set the timer for hours/minutes (H/M), push [MIN] again. Proceed as described above. The maximum measurable time is 99 hours and 59 minutes.
6. To reset the timer - Press [MIN] and [SEC] simultaneously to reset timer to "00:00".

Press [MIN] again to switch from timer to clock. The previously set timer will be recorded if it has not elapsed or been reset and will be shown again automatically when switched back to timer function.

brannan.co.uk

Manufacturer: S.Brannan & Sons, Leconfield Industrial Estate, Cleator Moor, Cumbria, UK, CA25 5QE
EU Authorised Representative: Comply Express Express Unipessoal Limitada, StartUp Madeira, EV141, Campus de Penteada, 9020 105 Funchal, Portugal



BRA.2022.REV1



Measurement for life.

Brannan